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EN I Heat vest with vibration massage function




INSTRUCTION MANUAL

Thank you for purchasing our Heat vest with vibration massage function. Please read this manual before use to ensure that you can enjoy your product for a long time and make full use of it.

Instructions

1. read all instructions.
2. Do not use during sleep.
3. Do not use on an infant.
4. Do not use this product on and by an invalid, sleeping or unconscious person or child. Persons with poor circulation must be carefully supervised during use.
5. switch off the product when not in use and unplug the power cord from the wall socket.
6. keep dogs, cats and other pets away from the heat waistcoat with massage function far away.
7. Use indoors only.
8. the device has a heated surface. People who are sensitive to heat must be careful when you use the machine.
DO NOT USE ON SENSITIVE AREAS OF SKIN.
9. burns can occur regardless of the control setting, please check your skin under the pillow frequently.
10. never use the cushion without the cover.
11. do not use clips/pins or other metallic means to fix the cushion.
12. do not sit on the cushion and do not crush it - avoid sharp folds
13. never lift the padding by the power cord alone.
14. do not use the cable as a handle.
15. check the inner cover carefully before each use. Do not use the padding if the inner cover shows signs of damage.
16. only loosely wrap the loop cord during storage, a tight wrap can damage the cord and the inner parts. This appliance must not be used by persons who are sensitive to heat and other persons who are not able to react to overheating.
17. The appliance must not be used by young children under 8 years of age unless they are supervised by a parent or guardian.
18. children under three years of age must never use this appliance, as they are not capable of reacting to overheating.
19. 19 When storing the appliance, allow it to cool down before folding it up.
20. 20. do not crumple the appliance or place objects on it during storage.
21. 21 Check the unit frequently for signs of wear or damage. If there are signs of damage or if the unit has been misused or does not work, return it to the supplier before further use.
22. This machine is not intended for medical use in hospitals.
23. never use when it is wet. **DANGER:** To reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions.
24. use only with the BLV-C09 controller.
25. use for too long can cause skin burns.
26. for prolonged use, should be set at a low level.
27. the switch or control unit must not get wet during washing and the cable must be positioned during drying so that no water runs into the switch or control unit.
28. if the power cord is damaged, it must be replaced by the manufacturer, the customer service or qualified persons to avoid any risk.
29. **CAUTION:** To avoid any risk of accidental resetting of the thermal cut-out, this appliance must not be powered by an external switchgear such as a timer or connected to a circuit that is regularly switched on and off by the utility company.
30. **If you do not use the waistcoat, please store it very carefully. Let the waistcoat cool down before folding it. Roll the waistcoat carefully or fold it loosely and without pressure to avoid damage to the heating elements.**
31. Do not place any objects on the waistcoat. Keep these instructions in a safe place.

INSTRUCTION MANUAL

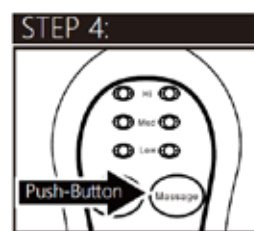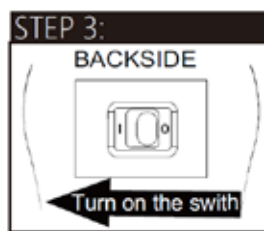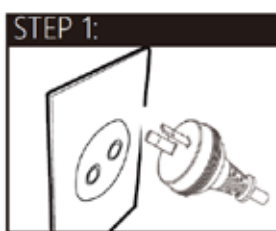STEP1: Insert the plug into a 220V-240V AC socket.

STEP2: Place the heat vest with vibration massage on your shoulders. Now connect the two front cushions with the protruding black velcro.


STEP3: Turn on the switch on the back of the controller

STEP4: then select the desired heat or vibrating massage function for the front of the controller.

IMPORTANT:

If the heat indicators light up or remain on, unplug the power cord from the wall socket immediately.


Preheat:

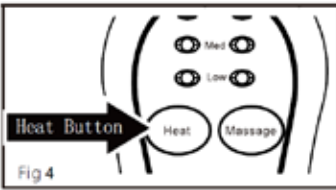Wait a few minutes for the heat waistcoat to preheat. **IMPORTANT:** If a red „HIGH“ lamp flashes, it means that an error has occurred. Immediately unplug the appliance and check the condition of the cable and the operation of the mains socket. If the flashing stops, stop using the waistcoat and contact customer service.


To add heat, turn on the switch on the back of the controller, then press the heat activation button and select the desired temperature: high or low.

Press 1 X for low heat
Press 2 X for medium heat
Press 3 X for lots of heat
Press again for the off function

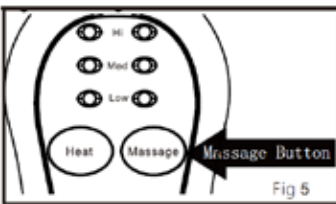To activate the vibration massage function, turn on the switch on the back of the controller, then press the vibration massage activation button and select your desired speed:

Press 1 X for slightly pulsating
Press 2 X for slightly permanent
Press 3 X for medium pulsing
Press 4 X for medium permanent
Press 5 X for strong pulsing
Press 6 X for strong permanent
Press again for the off function

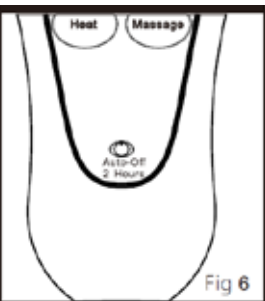
AUTOMATIC SWITCH-OFF

The auto-off switch is in the control unit on the Heat vest with vibration massage function. The Auto-Off function allows to save energy by automatically switching off after 2 hours of continuous use. The massage function switches off automatically after 15 minutes of use. **IMPORTANT:** The auto-off timer is reset after a power failure.

WASHING INSTRUCTIONS:

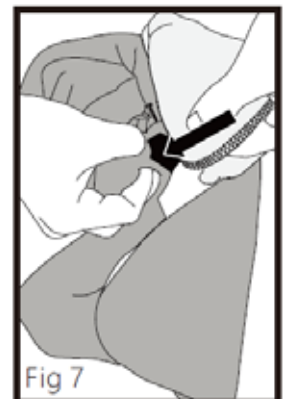ALWAYS DISCONNECT THE HEAT VEST WITH VIBRATION MASSAGE FUNCTION FROM THE MAINS BY PULLING THE PLUG! **NOTE: ONLY THE COVER IS WASHABLE! DO NOT WASH THE HEATING PAD! DO NOT DRY CLEAN! DO NOT USE BLEACH! DO NOT WRING OUT TO DRY.**

1. open the zip on the side and loosen the inner straps of the heating pad up to the cover. (See fig. 7). Then carefully remove the cover.
2. machine wash the cover in warm water on the delicate wash cycle and tumble dry on medium heat, then remove directly from the tumble dryer.
3. wash the cover by hand in warm soapy water. Simply hang it up to dry.
4. after drying, put the heating pad back into the cover and fasten it again.
5. close the straps with a zip fastener afterwards.

-  Do not bleach
-  Do not clean chemically or dry
-  Do not iron
-  Do not wash electric blanket and vibration element
-  only horizontal drying

Do not use metal needles

Do not use on children
under 3 years of age

Do not fold or crease

Read the instructions
carefully

ATTENTION: ALWAYS HANDLE PRODUCT WITH CARE. ALWAYS LAY FLAT OR FOLD CAREFULLY.
BENDING THE HEATING WIRES CAN CAUSE A DEFECT.
FOR DOMESTIC USE ONLY

100% polyester
66 x 85cm
220-240V~
50/60Hz 100Watt

CE

CE Declaration of Conformity

Mobiset GmbH hereby declares that this product is in compliance with the basic requirements of the applicable European directives. The declarations of CE conformity are available at www.mobiset.de.

General information

Copyright

This document is protected by copyright. Any duplication or reprinting, even in extracts, as well as the reproduction of the illustrations, even in modified form, is only permitted with the written consent of the manufacturer.

Notes on environmental protection

-  The packaging materials used are recyclable. Dispose of no longer needed Packaging materials according to local regulations.

-  This product must not be disposed of with normal household waste within the European Union will be. Dispose of the device at the municipal collection points.

Warranty/warranty

The warranty period is 24 months from the date of purchase. This includes a 6-month guarantee from the date of purchase.

The territorial scope of the warranty protection is Europe-wide.

The proof of purchase is valid as proof of this. The product has been thoroughly checked for quality and functionality before dispatch.

The guarantee is excluded in the following cases:

- in the event of modifications to the system without our approval
- in the case of use other than that stated above
- in the event of damage caused by falling

When making a claim under the guarantee/warranty, please return the complete product in transport-safe packaging to the point of sale from which the product was purchased. Describe the detected defect in the product and enclose the original proof of purchase with the shipment.

The seller reserves the right to remedy the defect before refunding the purchase price.

Necessary transport costs will be reimbursed to the customer if a justified claim under the statutory guarantee/warranty is established.

The use of the legal rights for defects is free of charge. These are not restricted by the guarantee.

Intended use

This device is only intended for use as a thermal waistcoat with massage function.

Any other or additional use is considered improper.

WARNING

Dangers can arise from the device if it is not used as intended and/or used in any other way. Use the device exclusively for its intended purpose.

Observe the procedures described in these operating instructions.

Claims of any kind for damage resulting from improper use and/or other use

Use is excluded.

The operator bears the sole risk.

Limitation of liability

All technical information, data and instructions for connection and operation contained in these operating instructions correspond to the latest status at the time of printing and are based on our previous experience and knowledge to the best of our knowledge.

No claims can be derived from the information, illustrations and descriptions in this manual.

The manufacturer accepts no liability for damage due to non-compliance with the instructions, improper use, improper repairs, unauthorised modifications or the use of unauthorised spare parts.

Safety

IMPORTANT SAFETY INSTRUCTIONS

PLEASE READ BEFORE FIRST USE

CAUTION: ALWAYS HANDLE PRODUCT WITH CARE. ALWAYS LAY FLAT OR FOLD.
SQUEEZING THE HEATING WIRES CAN LEAD TO MALFUNCTIONS.

IMPORTANT SAFETY INSTRUCTIONS

Please keep the operating instructions in a safe place for future use.

Warning: To avoid the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions.

This appliance can only be used by children under supervision.

Do not allow children to play with the appliance, cleaning and maintenance must be carried out by the user.

WARNING

The device must not be opened, thrown into fire, immersed in liquids or short-circuited.
will be.

Basic safety instructions

Check the unit for visible external damage before use. Do not operate a damaged or dropped unit.

Persons who are not able to operate the unit safely due to their physical, mental or motor skills may only use the unit under the supervision or instruction of a responsible person

Allow children to use the appliance only under supervision.

Have repairs to the appliance carried out only by authorised specialist companies or the customer service department. Improper repairs can cause considerable danger to the user. In addition, the warranty will become void.

Protect the unit from moisture and the penetration of liquids. Avoid direct sunlight.

Requirements for the installation site

For safe and fault-free operation of the unit, the installation site must meet the following requirements:

The unit must be set up on a firm, flat and horizontal surface or fixed to the wall using fixing materials.

Do not set up the device in a hot, wet or very damp environment or near inflammable material.

Before first use

Unpacking

Remove the machine and the operating instructions from the carton.

Remove all packaging materials and protective films from the unit before initial use.

NOTE

Keep the original packaging during the warranty or guarantee period of the device in order to send the device properly packed in case of warranty or guarantee claim. Damage in transit leads to the expiration of the warranty or guarantee claim.

Check scope of delivery

Check the contents for completeness and for visible damage.

Report an incomplete or damaged delivery to your supplier immediately.

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